

Cervical Spine Ankylosis Study

How to read the cases — a short picture guide.


There is **nothing to install**. The images stay on your own computer and **nothing is sent over the internet**. If you get stuck at any point, stop and contact us rather than guessing — nothing you do can damage the images.

1 Put the files together in one folder

You will receive one file called **Ankylosis-Study.html** and several files called **CASES-batch-1.zip**, **CASES-batch-2.zip** and so on.

1. Make a new folder on your Desktop called **ANKYLOSIS STUDY**.
2. Move **Ankylosis-Study.html** into it.
3. Right-click each zip file, choose **Extract All...**, and extract it into that same folder.

All the batches join together into one folder called **CASES**. When you are done it should look like this:

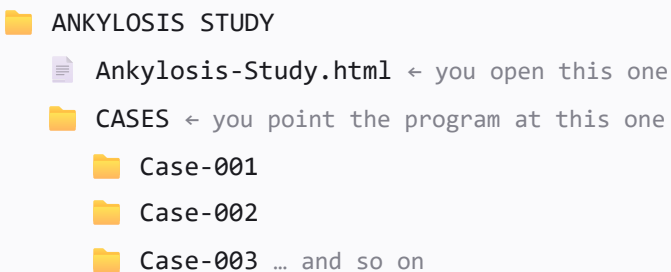

```
ANKYLOSIS STUDY
├── Ankylosis-Study.html ← you open this one
├── CASES ← you point the program at this one
│   ├── Case-001
│   ├── Case-002
│   └── Case-003 ... and so on
```

You can delete the zip files afterwards if you need the space.

2

Open the program

Double-click **Ankylosis-Study.html**. It opens in your web browser. Please use **Google Chrome** or **Microsoft Edge**. You will see this:

Cervical Spine Ankylosis Study

Click the button below to begin.

Open the image folder

Choose the folder named **CASES**. Your browser will ask whether to upload the folder — this is a standard message, nothing leaves your computer.

Click **Open the image folder**, then choose the folder called **CASES** and click **Select Folder**.

You will then see a message asking whether to upload the folder. This wording comes from the browser itself, not from us. Nothing is uploaded anywhere — it is only asking permission to read files already on your computer. Click **Upload** to continue.

The first case appears after a few seconds.

Case-001 Case 1 of 87 ————— 0 of 87 completed Save results file

Axial **Sagittal** Coronal **Bone** Soft tissue Reset view



ANKYLOSIS
Is ankylosis present at each level?

C1	Yes	No	Can't tell
C2	Yes	No	Can't tell
C3	Yes	No	Can't tell
C4	Yes	No	Can't tell
C5	Yes	No	Can't tell
C6	Yes	No	Can't tell
C7	Yes	No	Can't tell

Image < > Image 257 / 512

Brightness Contrast Zoom

Next case

- 1 **Top left** shows which case you are on and how many you have finished.
- 2 **Axial / Sagittal / Coronal** switch between the three planes. You can move freely between them.



Axial



Sagittal

- 3 **Bone / Soft tissue** set the standard window settings.
- 4 **The sliders under the picture** move through the images and change brightness, contrast and zoom.
- 5 **The right hand panel** is where you record your findings.

Moving around the images:

Move through the images

Drag the **Image** slider, or click   (hold to keep going). The mouse wheel and   keys also work.

Adjust brightness / contrast

Drag the **Brightness** and **Contrast** sliders. Or hold the left mouse button on the picture and move the mouse.

Zoom in and out

Drag the **Zoom** slider, or hold  and roll the mouse wheel

Slide the picture around

Hold the **right** mouse button on the picture and move the mouse

Put everything back to normal

Click **Reset view**, or **Bone** / **Soft tissue** for the standard settings

4 Recording your findings

For every level from **C1** to **C7**, choose one of three answers:

Yes

Ankylosis is present at this level.

No

Ankylosis is not present at this level.

Can't tell

This level cannot be judged — for example it is not fully included in the scan, or artefact obscures it. Please use this rather than guessing.

Then choose how confident you are about the case overall (**Low**, **Moderate** or **High**).

The  button only turns blue once every level has been answered. Anything still missing is listed underneath it. To change an answer, simply click a different button. To revisit an earlier case, click **Back**.

5 Saving and sending your work

Your answers are saved automatically as you go. You may close the page at any time and come back later — it will continue from where you stopped.

At the end of every session, click “Save results file” at the top right. A file is placed in your **Downloads** folder with a name like **ankylosis-FNGHKG-20260918-1430.csv**. Please email that file back to us.

It does not matter if you send it more than once, or if you have not finished all the cases yet. Just always send the most recent file.

? If something goes wrong

It says no images were found

The wrong folder was chosen. Click the button again and select the folder named **CASES** — not one of the **Case-001** folders inside it, and not a zip file.

The page is blank, or it opened in the wrong program

Right-click **Ankylosis-Study.html**, choose **Open with**, and select **Google Chrome** or **Microsoft Edge**.

It takes a while to open a case

Normal for the first case. After that, the next case is quietly prepared in the background while you are reading, so it should appear straight away.

I think I lost my answers

Look in your **Downloads** folder for a file whose name starts with **ankylosis-backup**. Open the program and click **Restore from a backup file** at the bottom of the first screen.

Anything else

Please contact us before carrying on, and mention the **six-letter code** shown at the top right of the screen.

Thank you for taking part. Please read the cases in the order the program presents them, and do not discuss individual cases with the other readers while the study is running.